

Panchakarma Products

▶ Taap Swedan Yantra (Dry Heat, Fomentation)

We have quality products in Heat Swedan Units. All our models are efficient and display quality in Panchakarma (Ayurvedic) products.



IMI - 2246: DELUXE SAUNA CABIN (for Two Persons)



Prefabricated, royal wooden cabin, size 60"x72"x75" high is suitable for 2 persons, thermally insulated ceiling & wall sections are made out of teak wood with optimum characteristic for use in high temperature and gives cosy atmosphere.

The sauna cabin is supplied with 6 Kw Stove & stove guard, Control Panel, Internal light, Rocks, Two layer Bench, Digital thermometer, Bucket & ladle. No extras are needed.

(Three Phase 440V.AC power supply is required, and Room will be assembled at destination.)

[\[Request Quotation / More info...\]](#) ↗

IMI - 2266: VIDYUT (I.R) BATH CABINET




Infra-Red bath cabinet, fitted with Eight I.R Lamps suitable for full body fomentation. Cabinet with thermostate to control the temperature is made out of good quality hardwood & ply, seat height is adjustable to need of the patient.

Wood polish finish.

[\[Request Quotation / More info...\]](#) ↗

IMI - 2559: INFRA-RED LAMP (Floor)

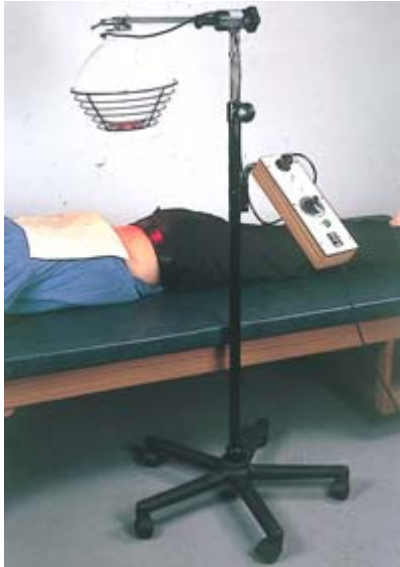

Floor stand model lamp for IR radiation treatments. IR bulb of 150w is fitted in metal casing with wire guard, automatic cut off timer of 15 minutes and an intensity control switch. Adjustable arm permits wide range of adjustments. Mounted on wide base with ball bearing casters.

[\[Request Quotation / More info...\]](#) 

IMI - 2565: ULTRAVIOLET LAMP (Floor)

No picture available

Floor stand model UV Lamp (Luminous) for Ultraviolet treatments. The UV bulb of 300w is fitted in metal casing with wire guard and automatic cut off Timer of 15 minutes. Adjustable arm permits wide range of adjustments. Mounted on heavy duty base fitted with ball bearing casters. Comes with suitable power cord.

[\[Request Quotation / More info...\]](#) 

IMI - 2573: SHORTWAVE MEDICAL DIATHERMY-500w. (with Fan)

High power medical diathermy incorporating a high power transmitting oscillator Valve & solid state rectifiers; provides R.F. power upto 500 watts for intense deep tissue heating with out discomfort at the surface. Housed in an elegant MOBILE CABINET & designed to meet heavy work load in a hospital.

SPECIFICATION:



- Frequency: 27.12MHz.
- Output: Continuous, maximum: 500 watts.
- Intensity: Adjustable in 5 steps.
- Electrodes: One pair Condenser pads with felt spacer & cover. Size: 15cm x 25cm.
- Tuning: Manual, with tuning control knob.
- Timer: 30 minutes timer with Auto-cut.
- Cooling: With 10cm dia. instrumental fan.
- Voltage: 220-240V; Selectable with inbuilt voltage selector switch.:

-Accessories:

- a) One pair condensor Pads cable.
- b) One power Main cable.
- c) Spare Fuses.
- d) One Shortwave Test Tube.
- e) Two Velcro Strap to tie / hold condensor pads with the body.

NOTE:

VALVE OF SWD-500 CARRIES NO GUARANTEE AS IT IS JUST LIKE AN ELECTRIC BULB THAT MAY FUSE ANY TIME. MUST USE VOLT. STABILISER WITH UNIT.

[\[Request Quotation / More info...\]](#) 